


Accessories – Transfer Stations

Step Turn Station

The step turn station is mounted outside a StoreX system and can rotate each microplate by an angle specified by the software. Since the StoreX handler is transporting the microplate along its longer axis, this feature may be necessary when an external handling system requires the microplate in another format.

The components used are the most reliable and of the highest quality available in the market. All manufactured components are electro-coated to protect against corrosion and contamination from external environmental conditions.

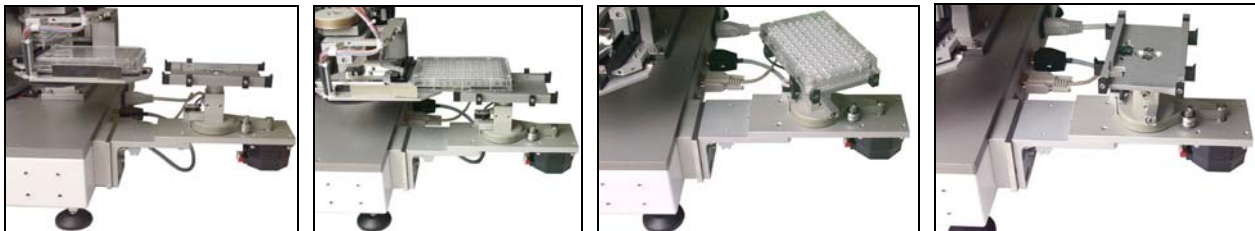

Order information

	<u>Order Nr.</u>
Step Turn Station	9111 11 21

Specifications

Plate Positions	1
Sensor	Plate presence sensor, Initialisation switch
Layout	Nest
Movement	Variable angles ranging from 0 degree to 300 degree
Application	Handling of plate provided at different angles